

TISSUE SELECTIVE THERAPY – TST33-S180

Designed for the treatment of Grade II–IV hemorrhoids, rectocele and rectal prolapse requiring partial or circumferential resection.

Precision Resection.
Enhanced Visibility.
Improved Control.

Mega-Windows™

Provides direct visualisation for precise tissue selection and targeted resection

Large Specimen Housing

Accommodates larger tissue volumes, expanding procedural versatility

ARL (Anti-Retrieving Lock)

Provides high-strength compression for secure tissue capture

PPC (Permanent Parallel Control) Technology


Ensures consistent staple formation and reliable haemostasis

Selective Sampling Accessories

Provide greater control over tissue selection for tailored resection

ASL Safety Lock

Helps prevent misfire and unintended re-firing for enhanced procedural safety

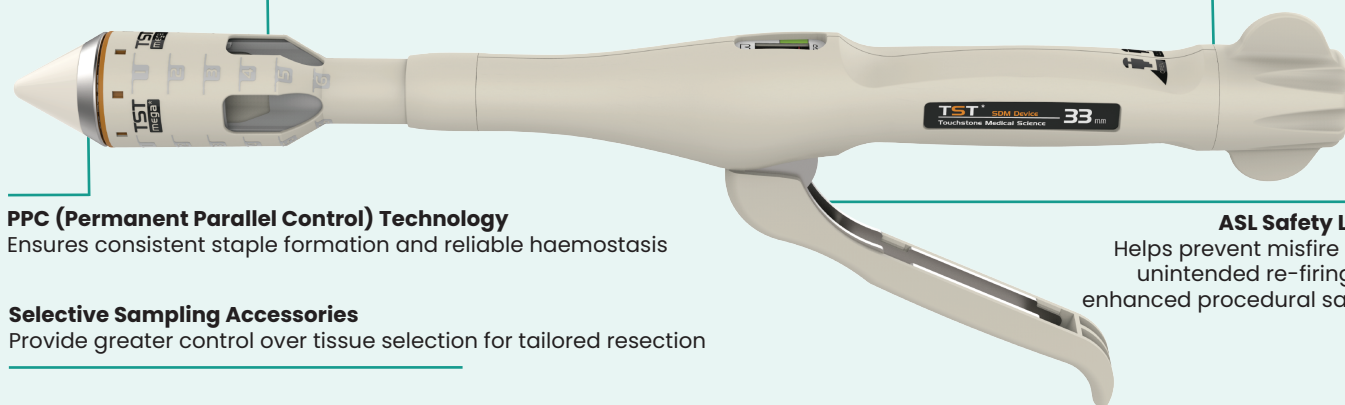

Stapler Code	Anvil Diameter (mm)	Staple Height (mm)	Staple Height Closed (mm)	Staple Housing Volume (ml)	Staple Quantity
TST33-S180	33	4.0	0.75–1.5	19	32
Accessory	Description				Pack Quantity
MEGA-DAO	Procto-Exam Kit with Anoscope, Obturator, Dilator & Dilator with Butterfly Wings 33mm				10

NEED MORE INFORMATION

Please contact us at
enquiries@lawmed.co.uk or **+44 (0)1932 260 838**.




Lawmed Ltd
Building 5 The Heights
Wellington Way
KT13 0NY
UK




www.lawmed.co.uk

TISSUE SELECTIVE THERAPY – TST33-S180

Segmented Purse-String

Circumferential Purse-String

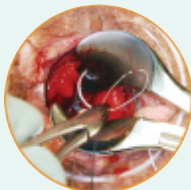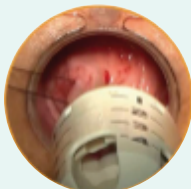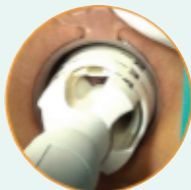- 1  Dilatation of Anus
- 2  Insertion & fixation of an appropriate Anoscope (Choose 1, 2 or 3 window)
- 3  Segmented Purse-String made
- 4  Split the suture to 2-3 strands and pull each strand into the traction hole separately with the suture threader
- 5  Turn the adjusting knob clockwise until the red indicator moves as far to the back as possible without over straining, then fire
- 6  Check the specimen

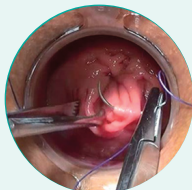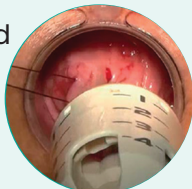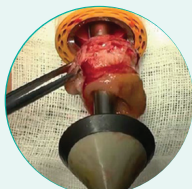- 1  Dilatation of Anus
- 2  Insertion & fixation of an appropriate Anoscope
- 3  Circumferential Purse String made
- 4  Insertion of TST-mega and tighten the suture on the centre rod
- 5  Turn the adjusting knob clockwise until the red indicator moves as far to the back as possible without over straining, then fire
- 6  Check the specimen

NEED MORE INFORMATION

Please contact us at
enquiries@lawmed.co.uk or **+44 (0)1932 260 838**.

Lawmed Ltd
Building 5 The Heights
Wellington Way
KT13 0NY
UK

www.lawmed.co.uk

TISSUE SELECTIVE THERAPY – TST36-S180

Designed for the treatment of Grade III–IV hemorrhoids, rectocele and rectal prolapse requiring larger volume partial or circumferential tissue resection.

Greater Capacity.
Single-Fire Efficiency.
Complete Control.

Mega-Windows™

Provides direct visualisation for precise tissue selection and targeted resection

35 cm³ Large Specimen Housing

Accommodates larger tissue volumes to support extensive prolapse resection

Barrier-Free™ Design

Provides additional space for larger tissue capture in a single firing

PPC (Permanent Parallel Control) Technology

Ensures consistent staple formation and reliable haemostasis

Selective Sampling Accessories



Provide greater control over tissue selection for tailored resection

ARL (Anti-Retrieving Lock)

Provides high-strength compression for secure tissue capture

ASL Safety Lock

Helps prevent misfire and unintended re-firing for enhanced procedural safety

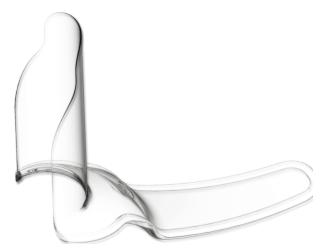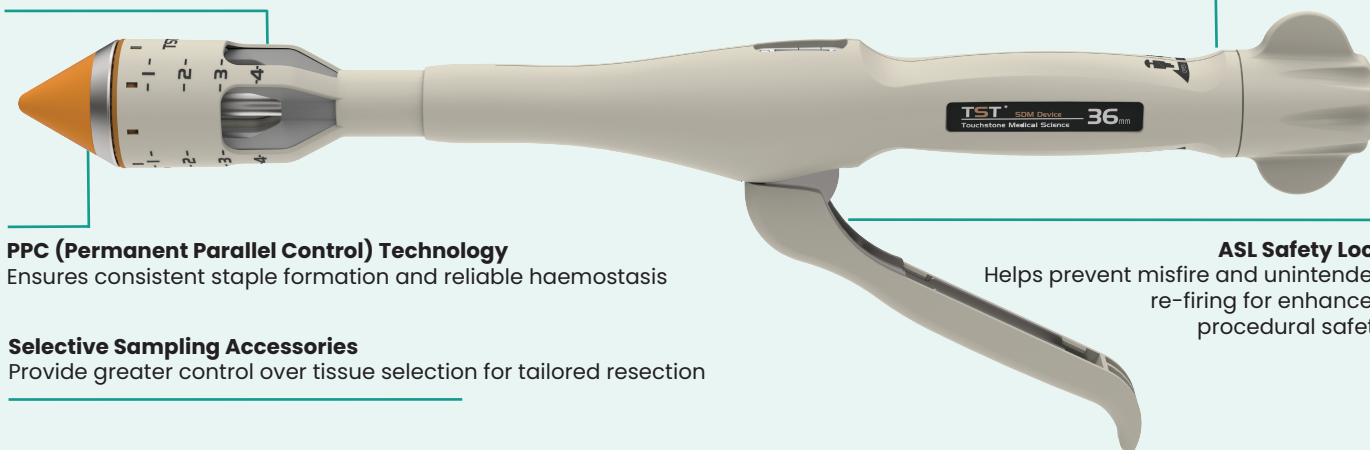

Stapler Code	Anvil Diameter (mm)	Staple Height (mm)	Staple Height Closed (mm)	Staple Housing Volume (ml)	Staple Quantity
TST36-S180	36	4.2	0.75–1.8	60	34

NEED MORE INFORMATION

Please contact us at
enquiries@lawmed.co.uk or **+44 (0)1932 260 838**.

Lawmed Ltd
Building 5 The Heights
Wellington Way
KT13 0NY
UK

www.lawmed.co.uk

TISSUE SELECTIVE THERAPY – TST36–S180


Parachute Purse-String

1




Anal inspection and dilation

3

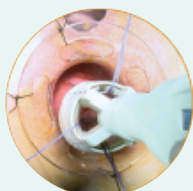

Split the purse string sutures into 4 strands. Pull each strand into the 4 traction holes separately with a suture threader. Adjust the depth of the purse string directly through the Mega-Windows

5

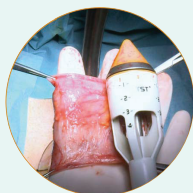


Examine the resected tissue and anastomosis

2



Make parachute purse string at 3, 5, 7, 11 o'clock

4

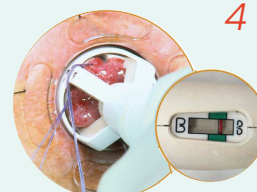

Turn the adjusting knob clockwise until the red indicator moves well into the green zone. The optimal position of the red indicator when in the green zone will be as far to the back as possible without over straining

NEED MORE INFORMATION

Please contact us at
enquiries@lawmed.co.uk or **+44 (0)1932 260 838**.

Lawmed Ltd
Building 5 The Heights
Wellington Way
KT13 0NY
UK

www.lawmed.co.uk