

References

- (1) Huo, Zhong-Chao et al. "Use of a disposable circumcision suture device versus conventional circumcision: a systematic review and meta-analysis." Asian journal of andrology vol. 19,3 (2017): 362-367
- (2) Cao, Y. J., et al. "Comparison of disposable circumcision suture device with disposable circumcision stapler and conventional circumcision." Chin J Clin (Electron Edition) 14 (2013): 6526-9.
- (3) Jin XD, Lu JJ, Liu WH, et al. Adult male circumcision with a circular stapler versus conventional circumcision: A prospective randomized clinical trial. Braz J Med Biol Res. 2015;48(6):577-582
- (4) Han, Hu et al. "Novel penile circumcision suturing devices versus the shang ring for adult male circumcision: a prospective study." International braz j urol : official journal of the Brazilian Society of Urology vol. 43,4 (2017): 736-745
- (5) Zhang, Guoyan et al. "Optimizing treatment strategies for pediatric phimosis and redundant prepuce: a comparative study of traditional circumcision and disposable circumcision stapler." Frontiers in pediatrics vol. 12 1394403. 22 Jul. 2024
- (6) Please refer to CCplus/SCC IFU

TOUCHSTONE



Advanced Stapler Circumcision Patient Information Guide

Distributed by



NEED MORE INFORMATION

Please contact us at
enquiries@lawmed.co.uk or +44 (0)1932 260 838
www.lawmed.co.uk

Advanced Circumcision Stapling

- A modern, shorter, stitchless procedures.
- less complications – faster return to work (1, 2, 3)



Shorter procedure time (vs. conventional circumcision)




Better cosmetic appearance, less complications



Faster return to work / daily activities

Comparison (1, 2, 3, 4, 5, 6)	Conventional Surgery	Stapler Circumcision
Duration	45-60 mins	10-15 minutes
Sutures/Stitches	Manual stitching required Stitch removal might be required	Silicone ring / staples Staples fall off spontaneously
Pain & Swelling	Moderate-high	Minimal to mild
Cosmetic Outcome	Dependent on manual cut	Clean, precise, symmetrical
Recovery Time	2-3 weeks	5-7 days

CCplus/Advanced Circumcision Stapling

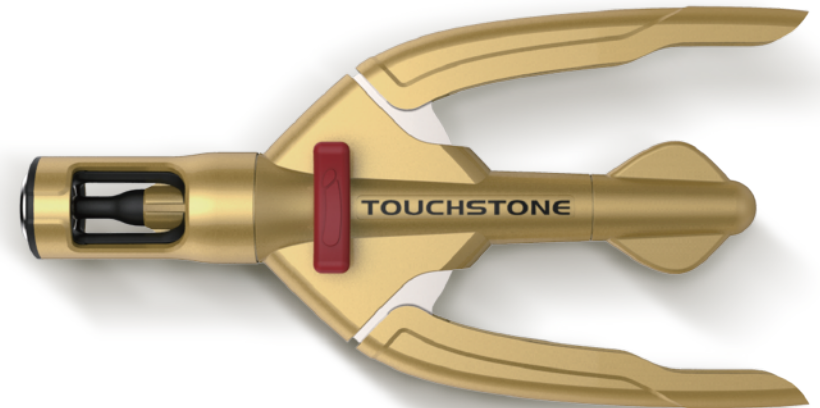


Mega Window

enhanced visual control,
precise positioning



Segregated Gasket

reduce the risks of bleeding and edema,
facilitates the healing process (4,5,6)

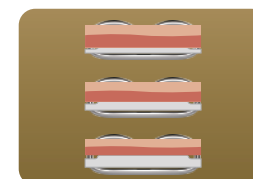

Optimized Staple Formation

B-shaped staples
reduce the risks of bleeding and edema,
eases the healing process (4, 5, 6)

Frequently Asked Questions Circumcision Stapling

Will there be pain during or after the procedure?

Local anaesthesia is applied prior to the procedure so you feel no pain during the process. After the local anaesthesia wears off you might experience some discomfort. Mild discomfort during recovery is manageable with prescribed over-the-counter pain relief.

What to expect after your circumcision?

Try to rest as much as possible during the first 3 days, keep the glans/penis facing upward. Try to maintain the penile area as dry and clean as possible, whilst the dressing is on.

Is bruising/hematoma/swelling normal?

Mild swelling, hematoma is normal after surgery and do not require any special treatment, resting can help to decrease swelling. If symptoms persist or become worse, please consult with your treating specialist.

When does the dressing come off and how do I care for my wound?

The dressing will be removed after 2-3 days (unless otherwise advised) by your doctor or nurse, but also can be removed at home. In case you remove it at home, moisten the dressing with a disinfecting solution (e.g. iodine) to facilitate the removal. After that, peel-off the dressing gently. Ensure to cleanse the wound area twice daily with an adequate disinfecting solution. Pad it dry gently with a sterile gauze. After dressing removal, if the tip is very sensitive, you can use some fragrant free moisturizer (e.g. petroleum jelly) for further protection.

When can I resume showering?

You can resume showering after 1 week. Gently washing the wound area can help the staples to fall off. Between washes it is important to keep the area clean and dry until healed.

When do the staples fall off?

The staples usually start to fall off by themselves after 2 weeks. Generally, it takes about 4 weeks to fall off completely.

When can I return to my normal activity?

After your surgery you can return to normal physical activities when you feel comfortable, unless otherwise instructed by your specialist. You may return to office work after approx. 4 days. However, avoid heavy or strenuous activity (including strenuous work) for the first couple of weeks. You can resume sport activities after 3 weeks, unless otherwise advised.

Avoid sexual activity (including masturbation) for at least 4 – 6 weeks as this may cause pain and bleeding.

Avoid spicy food and alcohol in the first weeks.



Key Facts



Pain will significantly ease about 2–3 hours after the anesthesia wears off. If you still experience strong discomfort, you may take oral pain medication.



Try to rest as much as possible during the first 3 days, keeping the glans penis facing upward.



Bandages can usually be removed 2–3 days after surgery. Before removing, moisten the gauze with iodine to gently peel it off. After removal, disinfect the area twice daily with iodine.



Mild swelling or hematoma is normal and does not require special treatment. Resting more can help reduce swelling.



You may shower one week after surgery. Gently washing the area can help the surgical staples fall off more easily. After bathing, make sure to keep the surgical site dry.



The staples usually start to fall off by themselves after 2 weeks, generally it takes about 4 weeks to fall off completely.



Avoid spicy foods and alcohol; ensure rest and recovery.



Avoid sexual activity (including masturbation) for at least 4–6 weeks.

Post-Circumcision Care Underwear For Children

Product info



Size

Weight Range

130

<25kg

140

<31kg

150

<36kg

Material

Main fabric: 95% cotton + 5% spandex

(Soft, breathable, moisture-wicking)

Protective cover material: TPU

(Strong support, durable, non-deforming)

Features Benefits

- TPU protection prevents wounds contact and reduces friction
- Magnetic opening allows observation and easy care
- Ventilation holes keep the area dry

Instructions

1. Wear like underwear
2. Adjust the protective pouch to cover the surgical area
3. Secure for all-day comfort