


Steps for Use

Powered Universal Handle*

1. Select the stapler size appropriate for the target surgical procedure and site.

Check the following;


- Package integrity
- Expiry date
- Sterilisation symbol

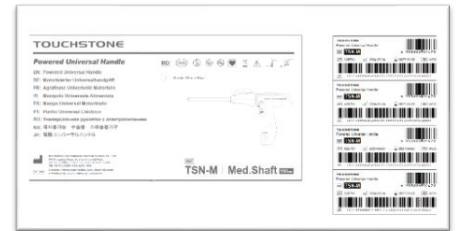

2. Using sterile technique, carefully remove the Instrument, Battery Pack, and Reload from their respective packages. Avoid flipping them into the sterile field to prevent damage.



3. Install the Battery Pack into the Instrument until it clicks into place. The Instrument will perform a self-check, and the battery indicator light will turn green, confirming proper installation and a full charge.

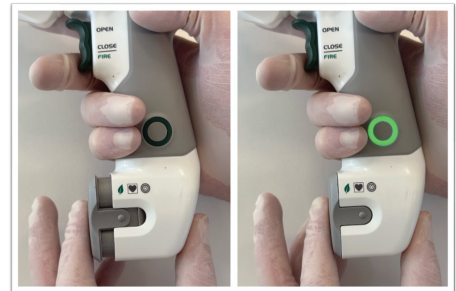

4. Attach the proximal end of the Reload to the Instrument shaft pin, ensuring that the alignment indicators on both components are matched (arrow to arrow).



5. Press the Reload onto the Instrument shaft, rotate it clockwise until it clicks, and wait for an audible beep. Then remove the orange shipping wedge.



Lawmed Ltd
Building 5 The Heights
Wellington Way
KT13 0NY
UK

www.lawmed.co.uk
enquiries@lawmed.co.uk
+44 (0)1932 260 838

Steps for Use

Powered Universal Handle*

6. The green status indicator light confirms successful installation and displays the staple line length.




7. Close the Reload jaw by pressing the Close/Fire Button down. An audible beep indicates that the jaw is fully closed.



8. Press and hold the Open Button until the jaw is completely open. An audible beep indicates that the jaw is fully open.



9. Close the Reload jaw by pressing the Close/Fire Button to facilitate insertion of the stapler into the patient cavity through a trocar sleeve, until the articulation joint is visible.

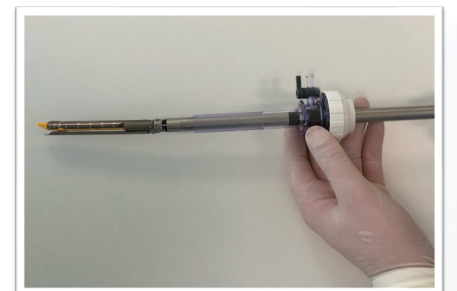

10. Open the jaw by pressing and holding the Open Button until the jaw is completely open. Position the target tissue between the Reload jaws, then clamp the tissue by pressing the Close/Fire Button down. Ensure the tissue is correctly positioned before fully closing the stapler jaw.



Lawmed Ltd
Building 5 The Heights
Wellington Way
KT13 0NY
UK

www.lawmed.co.uk
enquiries@lawmed.co.uk
+44 (0)1932 260 838

Steps for Use

Powered Universal Handle*

11. Press the green Safety Button. An audible beep will sound, and the status indicator light at the back of the Instrument will flash green.




12. To fire the stapler, press the Close/Fire Button down until the knife reaches the distal end. An audible beep indicates that firing is complete.



13. To open the stapler jaw, press and release the Open Button until the knife returns to its original position. Then press and hold the Open Button again until the jaws fully open and an audible beep is heard.



14. To uninstall the cartridge, insert the opening tip of the orange shipping wedge into the Reload anvil and slide it to the end. Press the tip of the anvil against the cartridge tip and push upward (toward the anvil) to release the cartridge from the Reload jaw.

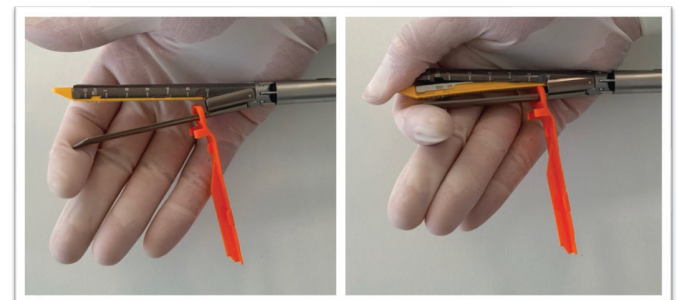

NOTE: Before reloading a new cartridge, hold the Instrument upright with the anvil and cartridge jaw fully submerged in a sterile solution. Swish briskly to remove any loose staples.



Lawmed Ltd
Building 5 The Heights
Wellington Way
KT13 0NY
UK


www.lawmed.co.uk
enquiries@lawmed.co.uk
+44 (0)1932 260 838

Steps for Use

Powered Universal Handle*

15. To load a new cartridge:

- Insert the cartridge of your choice at an approximate 30° angle toward the proximal side of the jaw until fully seated.
- Snap the cartridge down until a “click” is heard.

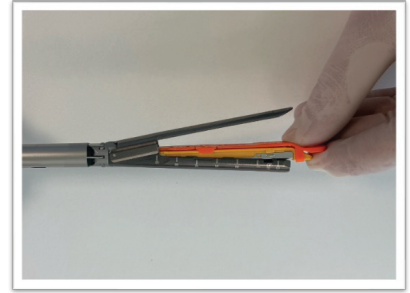

16. To unload the Reload, press the Unload/Unlock Button, rotate the Reload counterclockwise approximately 45°, and remove it from the Instrument.

To reload a new Reload into the stapler, repeat the Reload installation steps.




Caution: If pressing and holding the Open Button does not return the cutting knife to its initial position and the jaw remains closed, follow steps 17–19.


17. Ensure that the Battery Pack is securely installed and that the Instrument is powered on. Continue pressing and holding the Open Button. If the cutting knife does not return and the jaw does not open, use the Manual Override System.



18. To use the Manual Override System, remove the access panel labeled “Manual Override” on the top of the Instrument handle. The emergency lever will then be exposed.

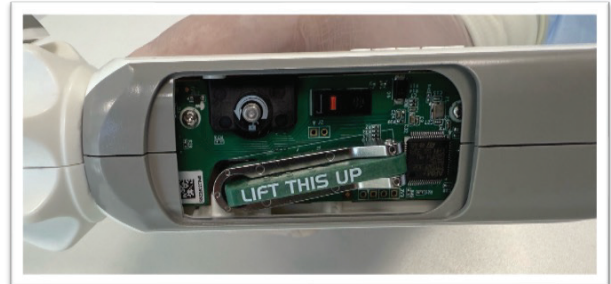

19. Lift the emergency lever, place it on the screw, and rotate it counterclockwise until the Cutting Knife returns to its initial position and the jaw is fully opened. Discard the Instrument.

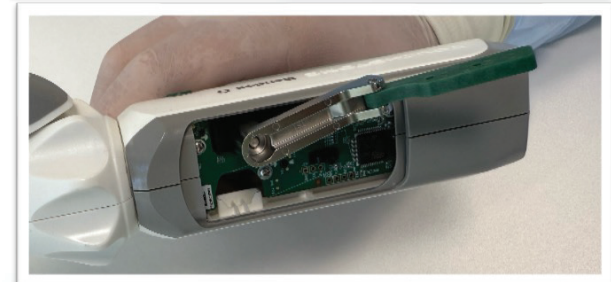

Lawmed Ltd
Building 5 The Heights
Wellington Way
KT13 0NY
UK

www.lawmed.co.uk
enquiries@lawmed.co.uk
+44 (0)1932 260 838