


**LINEAR CUTTER AND RELOADS**

Stable tissue positioning  
 Balanced closure  
 Secure staple lines

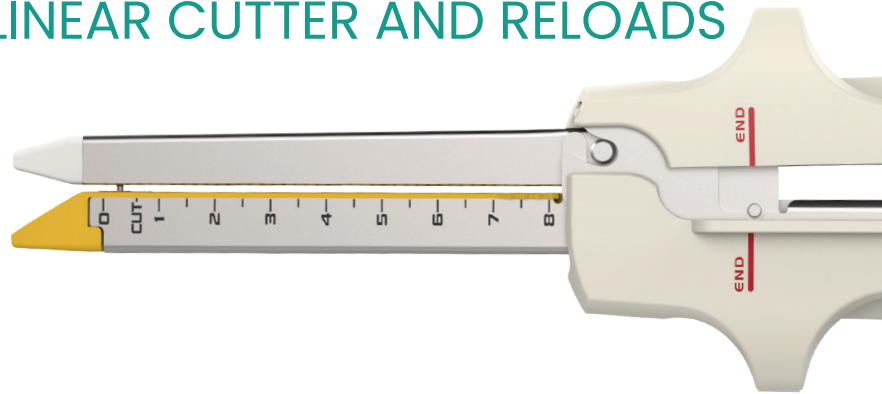

**BLC™ - Balanced Leverage Control System**

**WSL™  
 Wider Safety Lip**

Provides a wider surface to help


**Thickened Steel Construction**  
 Enhances instrument strength and stability

**Double Balancing Points**  
 Supports uniform tissue compression along the staple line

**Tissue Positioning Point**  
 Helps prevent tissue extrusion and consistent staple formation

**Tissue Positioning Chain**  
 Helps secure tissue within the instrument forks and minimizes tissue slippage

**Shielding Blade**  
 Covers blade before firing

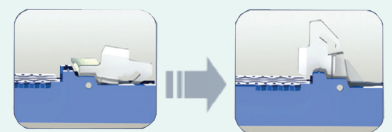

Before Firing

After Firing

**TPM™ - Tissue Positioning Management**

**NEED MORE INFORMATION**

Please contact us at  
**enquiries@lawmed.co.uk** or **+44 (0)1932 260 838**.

Lawmed Ltd  
 Building 5 The Heights  
 Wellington Way  
 KT13 0NY  
 UK

**www.lawmed.co.uk**

## LINEAR CUTTER AND RELOADS

**NEW KNIFE WITH EVERY FIRING** - Delivers consistent cutting performance with every reload

**FLOATING CARTRIDGE** - Adapts to tissue variations for more consistent compression

**UP TO 12 TOTAL FIRINGS** - One instrument. Up to 11 reloads, 12 total firings


**DUAL-SIDE LENGTH INDICATORS** - For precise positioning

**END-LINE INDICATION** - For accurate placement

**DOUBLE-SIDED FIRING LEVER** - Easy activation from either side

**FAST-RELEASE BUTTON** - Quick and easy instrument opening

**ANTI-SLIP FINGER PAD** - Secure grip and enhanced control

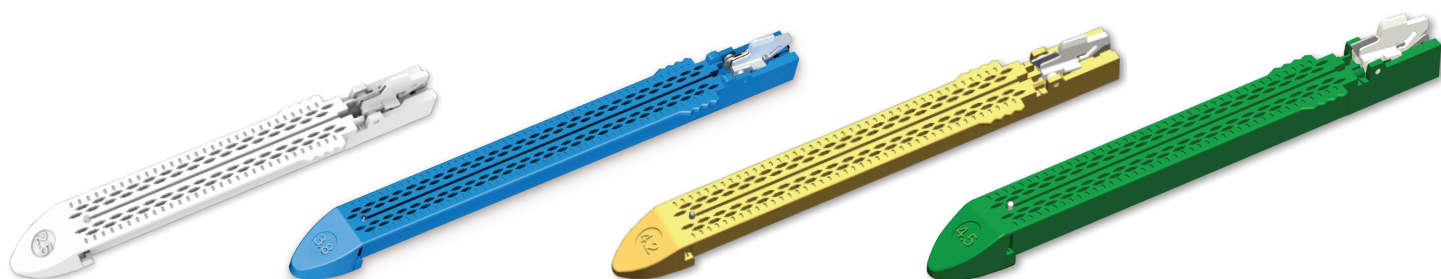

Stapler Code	Reload Code	Stapler Length (mm)	Staple Height (mm)	Staple Height Closed (mm)	Staple Quantity	Colour
LC6025	LCC6025	60	2.5	1.0	64	○ White
LC6038	LCC6038	60	3.8	1.5	64	● Blue
LC8038	LCC8038	80	3.8	1.5	84	● Blue
LC10038	LCC10038	100	3.8	1.5	104	● Blue
LC8042	LCC8042	80	4.2	1.8	84	● Gold
LC6045	LCC6045	60	4.5	2.0	64	● Green
LC8045	LCC8045	80	4.5	2.0	84	● Green
LC10045	LCC10045	100	4.5	2.0	104	● Green

### NEED MORE INFORMATION

Please contact us at  
[enquiries@lawmed.co.uk](mailto:enquiries@lawmed.co.uk) or **+44 (0)1932 260 838**.

Lawmed Ltd  
Building 5 The Heights  
Wellington Way  
KT13 0NY  
UK

[www.lawmed.co.uk](http://www.lawmed.co.uk)